

Stay safe when dealing with spills of Cytotoxic Drugs

Cytotoxic drugs are handled in areas such as pharmacy, wards, outpatients and even in the home. These substances are also transported between the various locations. A spill of a cytotoxic drug can therefore have serious consequences.

In the event of a spill a system is needed that completely manages the incident, keeping you safe and at the same time effectively removing the hazard and decontaminating the residue. Furthermore, the incident should be managed without causing alarm to patients.

Our spill system for Cytotoxic Drugs ensures a rapid and complete response in a low-key manner. The Spill Kit can be used as a permanent single location system or in a portable manner during transport.

Protective Equipment

2 full PPE sets includes face piece, eye shields, 2 pairs of gloves, lab coat and overshoes. All items selected for ease of use and 'low key' appearance to avoid causing alarm to patients.



Liquid clean up

T-Mop absorbent head rapidly absorbs and locks the liquid ready for waste disposal.

Protective absorbent kneeling mat and tweezers for glass also provided.



Highly visible

Forms an evacuation point containing everything needed to manage a spill and ensure a quick response.

Low profile design

Saves space in cramped areas. Manufactured from fire-resistant, anti-fungal, anti-bacterial wipe clean material.

Information where it is needed

The top compartment contains:

- ◆ spill instructions and full procedure
- ◆ warning signs
- ◆ report forms



Easy access sections

Allows **QuickBags** to be removed quickly and taken to the spill.



Contains 3

One Shot QuickBags

Each containing all equipment needed to manage a spill. QuickBags are available individually as a portable option for use in transport such as on trolleys or in vehicles.


Powder clean up

Sterile water to mobilise powder which is then absorbed by pads.



Waste disposal

Heavy-duty waste disposal bag with return labels makes sure that the waste is removed through the correct routes.

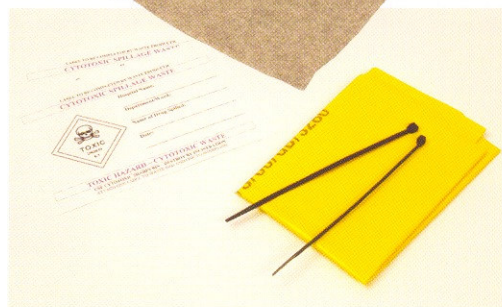

3-stage decontamination

1 pod of sterile water is used with 3 lint free wipes and this is repeated 3 times to ensure complete removal of residue.

Cairn Technology Ltd
Sheffield Technology Park
60 Shirland Lane
Sheffield, S9 3SP

Tel. 0114 244 9551
Fax. 0114 296 1005
info@cairntechnology.com
www.cairntechnology.com



Contents of the Cytotoxic Drug Spill Kit



Base Unit

Product code CT0500

Wall-mounted base unit containing 3 one-shot **QuickBags**



QuickBag

Product code CT0520

Containing all equipment to deal with a spill.

Contents of One Shot QuickBag

2 PPE sets containing face piece, eyeshields, 2 pairs of gloves, labcoat and overshoes

| Liquid clean up containing T-Mop, kneeling pad and tweezers

| Powder clean up containing pads and sterile water

| Decontamination containing sterile water and lint-free wipes

| Waste disposal containing waste bag, ties and returns labels

| Warning sign

| Report and management sheets

Our **Spill Systems** are now the choice of over 300 hospitals in the UK.

To place an order or request an information pack contact us now on:

Phone: +44 0114 244 9551 Fax: +44 0114 296 1005 E-mail: info@cairntechnology.com

About Cairn Technology

We provide consulting services and products that dramatically enhance personal and organisational performance where chemical and biological materials are handled. By teaming up with major NHS Trusts and Universities in the UK we have applied our technical knowledge to the daily challenges of handling chemicals, biological fluids and drugs.

Many of the NHS Trusts we have worked with have implemented our range of spill stations on a Trust wide basis ensuring a coherent and cost effective response for all sites and applications.

Visit us at www.cairntechnology.com for free downloads and to view our full range of products and services.

Cairn Technology Ltd
Sheffield Technology Park
60 Shirland Lane
Sheffield, S9 3SP

Tel. 0114 244 9551
Fax. 0114 296 1005
info@cairntechnology.com
www.cairntechnology.com

